


# Mobile Companion iMeasure.ai

## How-To Guide

### Logging into the App

- Upon entering the app, click the **Login** button
- When prompted, click **Continue** to access the login
- Enter your username and password for IntellicureEHR, and click **Sign In**

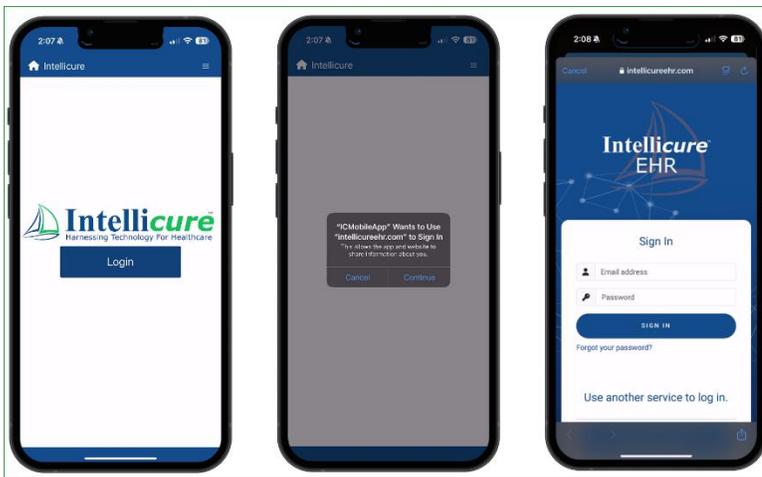


### Home Page

- The home page shows the user logged in and clinic
- Quickly access Photos or Problems using the menu bars



### Navigation

- Along the bottom of your screen is the **Navigational Menu** - Home, Photos, Patients, and Problems.

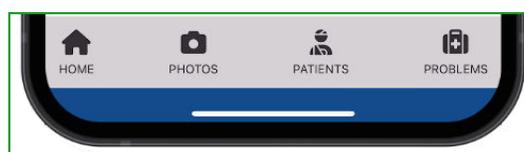

- In the top right corner of your screen, locate the menu icon that will appear on all tabs and pages after the home page

- **Options include**

- Home
- QR Search
- Change Databases
- Change Default Clinic
- Log Out



- The Mobile Companion supports several workflows. From Photos and Problems you have the ability to take a photo, enter measurements, document appearance, treatment orders, and add a new problem.
- The Patients tab is a future feature.


## Photos

- Navigate to **Photos** by clicking the **blue box** on the home screen, or **Photos** found in your navigational menu at the bottom of your screen




## Photos – Patient Search

- Here you can see a list of patients you are a clinical participant on, or search for any patient using Last Name, First Name, or Account Number.
- In the top right, turn the filter off to expand the list below.
- The **black L box** indicates the encounter is currently **locked**
  - You may have the ability to upload to a locked encounter (setting). You **cannot** add a problem to a locked encounter.

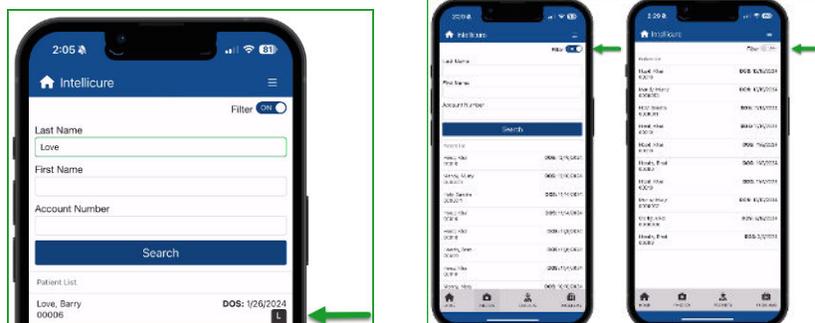


## Photos – Existing Problems

- Locate your patient and open the chart
- The patients name, date of birth, date of service, visit type, place of service, MRN, account number, and site are located at the top.
- Below you see a list of existing problems, click on the **Problem** to see further details
  - This section is called the **Photo List**

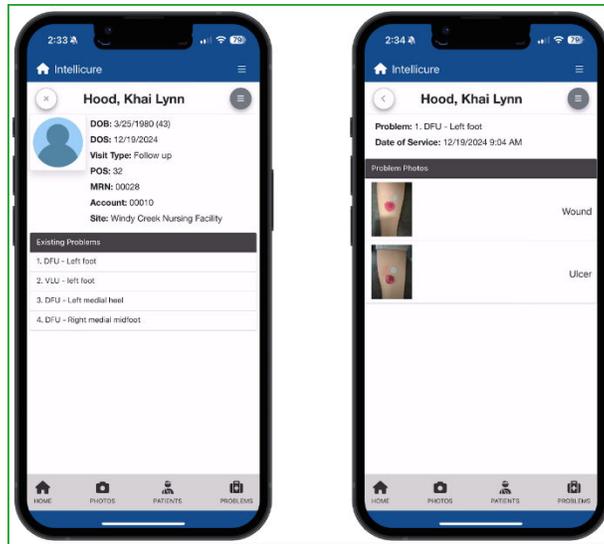


## Photos – Add a Photo

- From the **Photo List**, click the **menu icon** in the top right corner to add a photo to that problem

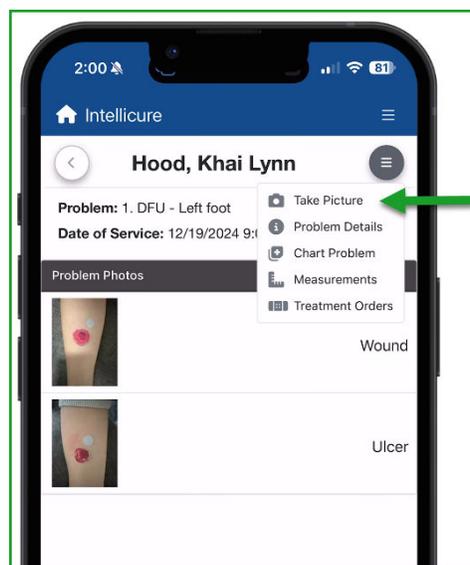


## Photos – Retaking a Photo

- Choose the photo that needs retaken, click the **menu icon**, choose **Retake Photo**
- After capturing, label the photo, and click **Save**

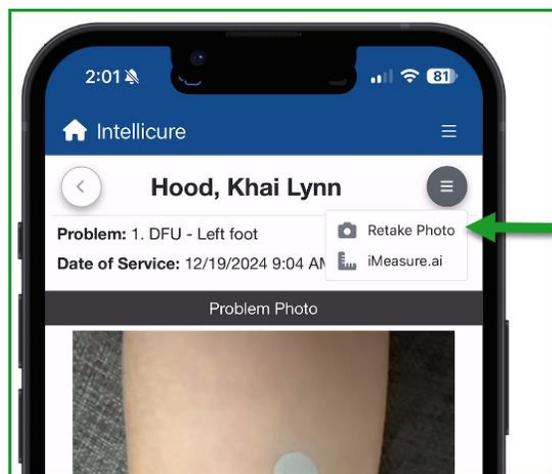


## Photos – iMeasure.ai

- To measure a wound, click the **Problem** in which the wound applies to in the **Existing Problems List**
- Click the **menu icon** in the top right corner to **Take Photo**
- Place your **gray marker** at the top of the wound about **1cm** from the wound edge using the cephalad to caudad or head to toe orientation.
  - This will allow you to consistently capture length and width problem to problem.
  - Make sure the sticker is not inside the wound margins, but right outside
  - Recommendation: lock screen rotation
- **Take Photo**, then **Label** the Photo and **Save**
- Once the photo is saved, **iMeasure.ai** will automatically measure the length and width of the wound.
- Use the **green points** to adjust the edges, if needed, then click **Looks Good**
- Select the deepest point of the wound (shown by **yellow point**), then click **Save Depth**

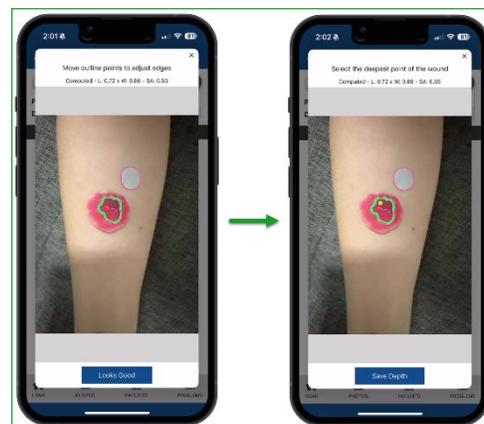


- Enter the wound depth and make any adjustments here, then **Save**
- To add additional measurements, click the **menu icon** and choose **Measurements**
  - Pre/Post Debridement
  - Sinus/Tunnels
  - Undermining

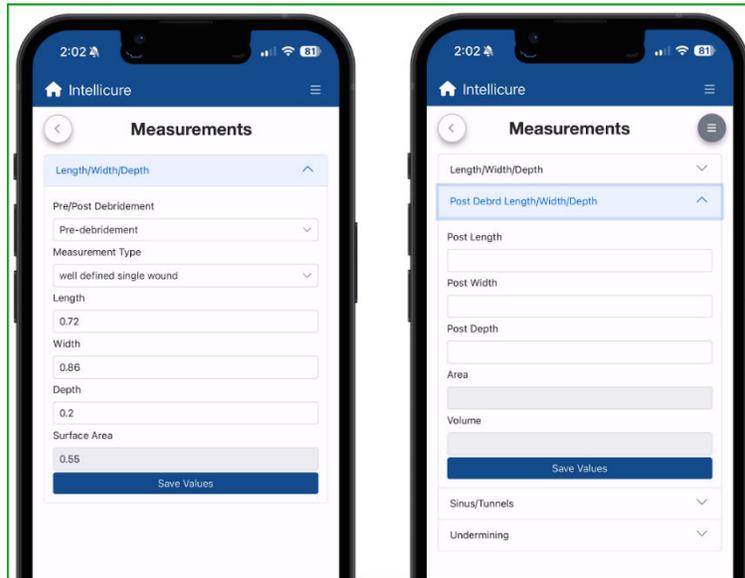


## Photos – Charting a Problem

- To chart appearance on the selected wound, click **Chart Problem** from the **menu icon**
- **Last values** will pull forward documentation from previous visits if selected
- Save documentation in each category before moving to the next.

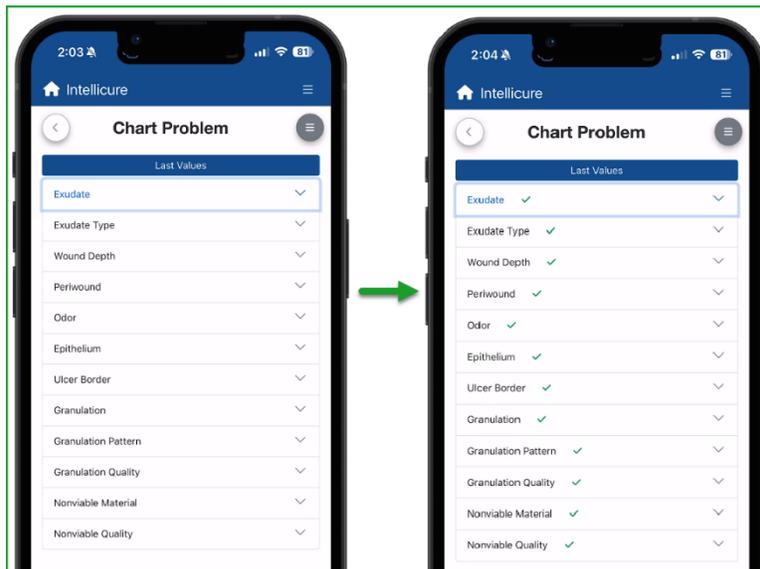


## Photos – Treatment Orders

- To document treatment orders from the Photos page, click the **menu icon** and choose **Treatment Orders**
- **Last values** will pull forward documentation from previous visits if selected
- Save documentation in each category before moving to the next.

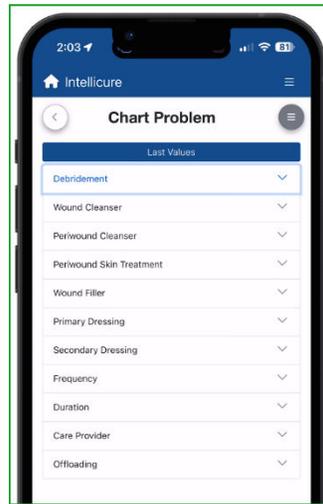


## Problems

- Click the **Problems** tab in the navigational menu and locate your patient
- Click the **menu icon** in the top right corner to add a new problem
  - Turn common problems on to see your customized list **OR** choose from the dropdown below

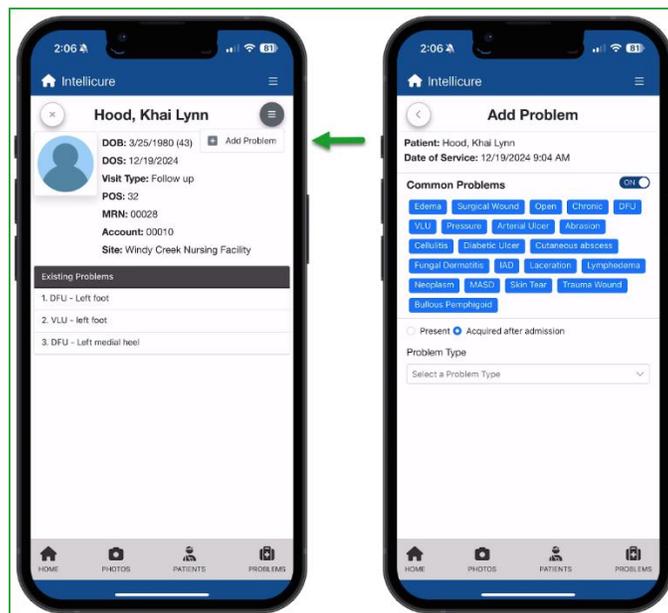

- Choose clinical terms that further describe the Problem
- As you choose terms, more will populate to ensure a detailed description
- The detailed description will populate the ICD-10 code
- Enter the **Location** and **Onset**, then **Save**
- The problem is now added to your existing problems list

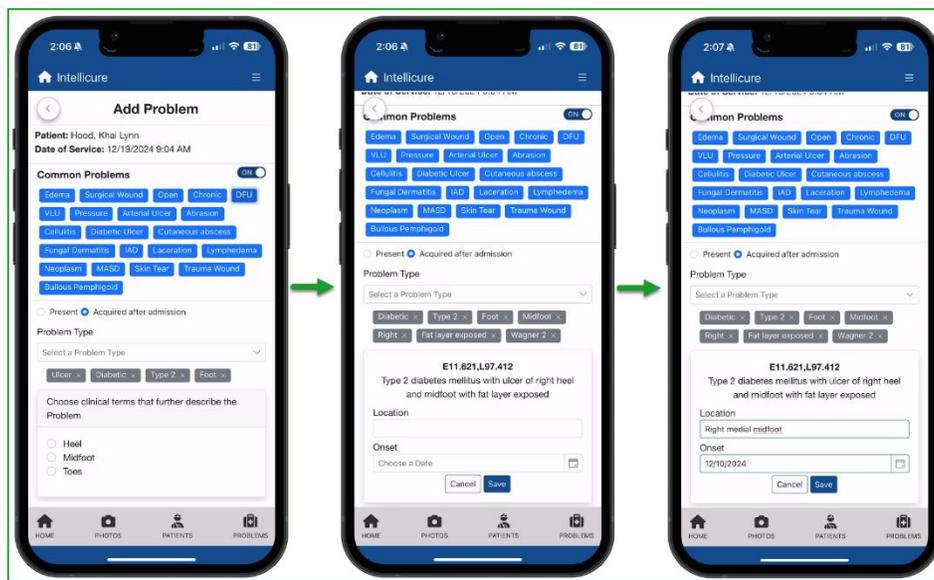

- Click on the **Problem** to see:
  - Problem Details
  - Onset Date
  - Noted Date
  - Diagnosis Codes
  - POA Status
- Click the **menu icon** in the top right corner to document:
  - Chart problem (appearance)
  - Measurements
  - Treatment Orders
  - Problem Photos
    - Take Photo

